

Cybersonic® Classic - The Original
Our most popular & successful model.

- World's first and only "harmonic resonance" cleaning technology
- Exclusive U.S. Patented sonic system
- Disintegrates bacterial plaque on contact
- FREE Replacement Brushes For Life*

*small S&H charge applies

One Year Limited Warranty*

(*FREE Lifetime Warranty with any active brush replacement program)

The Cybersonic® Classic Limited Warranty covers mechanical and/or electrical failures of the power system. Under warranty terms and conditions, the system shall be repaired or replaced at the manufacturer's discretion. Defective system and proof of purchase must be submitted along with warranty claims. This warranty does not cover:

- Brush heads
- Floss cartridges
- Damage caused by misuse, abuse, neglect or alterations
- Normal wear and tear, including chips, scratches, abrasions, discolorations or fading

Warranty claims and proof of purchase should be forwarded to:

AMDEN CORPORATION
P.O. Box 2940
Mission Viejo, CA 92690-0940
Toll-free number: 800-770-7525
www.sonictoothbrush.com

© 2008 Amden Corporation.
Cybersonic® and Sonic White® are registered trademarks of Amden Corporation.

40515

Classic
cybersonic®
high tech simplicity



**THE
ORIGINAL**

You asked for it and now,
back by popular demand!



cybersonic *Classic*®

Thank you for choosing Cybersonic® Classic, our most popular and successful product ever. It is The Original and world's only "harmonic resonance" sonic oral-care system. You get the same sonic technology, reliable performance, and ease-of-use that made Cybersonic® famous, plus with Cybersonic® Classic, we've added a few modern improvements to make it better than ever. Your system comes with The Original sonic toothbrush head, a new larger brush head, an improved sonic flosser with disposable cartridges, and a re-designed ergonomically shaped tongue cleaner for better cleaning.

If you missed our commercials, our exclusive U.S. Patented system works something like a musical tuning fork. In other words, the brush is specially tuned to "resonate" back & forth at sonic speeds, literally disintegrating plaque on contact!

Also, take advantage of another Classic favorite, which is our FREE Brush Replacement Program. Pay a nominal S&H fee and we'll automatically provide FREE brush heads every 3 months AND protect your Power Handle & Charging Stand with a lifetime warranty!


As always, we're ready to assist you or answer any questions via e-mail at: help@amdencorp.com, or you can call us toll-free at: 800-770-7525.

To your health,

Ronald K. Murayama, D.D.S.
Chairman & Founder



Your Cybersonic® Classic system includes:

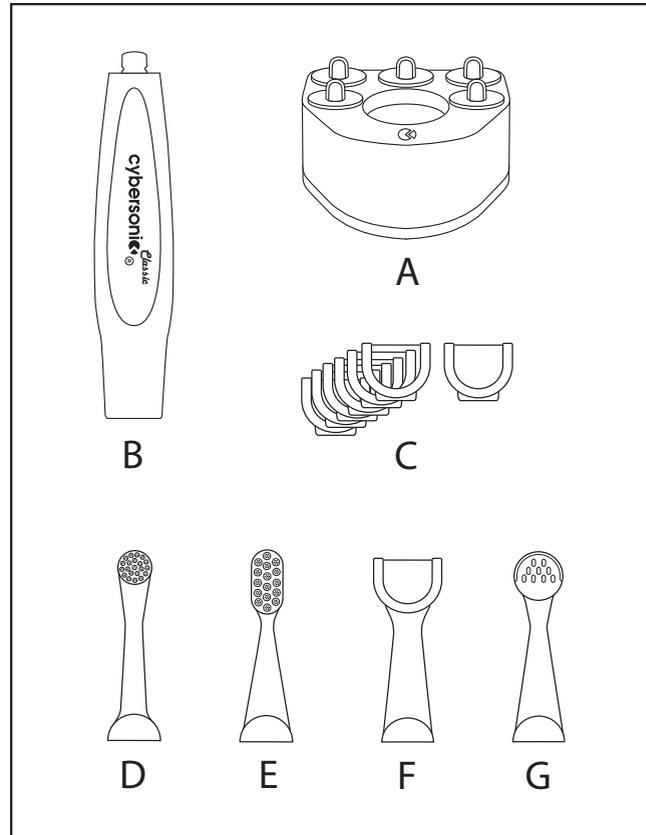

- A - Combo Charger & Accessory Stand
- B - Sonic Power Handle
- C - Disposable Flossing Cartridges
- D - Original Cybersonic® Sonic Brush Head
- E - Large Sonic Brush Head
- F - Sonic Flosser Attachment/Floss Holder
- G - Improved Sonic Tongue Cleaner

U.S. Patented Harmonic Resonance Technology

- Exclusive dentist-developed U.S. Patented sonic technology
- World's first and only "harmonic resonance" sonic system
- Disintegrates bacterial dental plaque on contact

Sonic Brush Heads

- 2 sizes accommodate the largest mouth or small kids
- Safe end-rounded Dupont® Nylon bristles
- Extra-soft and gentle .006" diameter brush bristles
- Brush head cap keeps bristles clean and sanitary

Power Handle

- Ergonomic design makes it easy to hold and control
- Hermetically sealed and 100% waterproof
- Can last one week on a single charge

Sonic Flosser

- Makes it easy to floss every day
- Pre-strung disposable floss cartridges
- Snaps on and off with ease
- Holds floss in the proper position

Tongue Cleaner

- Removes bad breath at the source, which is the tongue
- Works like a mini carpet cleaning agitator
- Low profile head and ergo-angle reduces gag reflex
- Freshens breath naturally instead of minty cover-ups

Charging Stand

- Small footprint uses minimal counter space
- Self-drains and easy to clean
- Adjust cord length on bottom of base
- Holds up to 6 attachments (including handle)

Sonic White® One-Step Bleaching Gel (OPTIONAL)

- Revolutionary one-step bleaching system
- Brushes on just like toothpaste.
- Kinetic energy stimulates the release of oxygen — the active whitening ingredient
- The easiest way to a brighter, effervescent smile
- Ideal for maintaining professional whitening treatments at home
- Active ingredient 7.5% Hydrogen Dioxide



Note: One syringe contains a one-week treatment; six days on, one day off. A one-week treatment is enough for normal yellowing or stains. Sonic White is safe for daily use. To order, call customers service toll free at 800-770-7525 or you can also place your order online at www.sonictoothbrush.com

Free Replacement Brushes & Lifetime Warranty

- Automatically receive FREE brush heads every 3 months
- Any active replacement program automatically warranties your Cybersonic® Classic system for life
- Just pay a small S&H
- Program is good for as long as you own your Cybersonic®
- Cybersonic® is the last toothbrush you'll ever need to buy!
- To sign up or for more information: www.sonictoothbrush.com or call toll-free at 800-770-7525

Cybersonic® Classic is as easy to use as it is to maintain. Please follow these simple instructions to get started and to ensure that your system works properly and effectively for years to come. If this is your first time using a power toothbrush, you'll notice a little tingling or "tickle" sensation. This is normal and minimizes after continued use.


1. Charge the Power Handle the first 24 Hours.

Before using your system for the first time, it is very important to charge the power handle for at least 24 hours before the first use. This also helps prolong the battery system.

- Plug charger into outlet
- Charging stand and power handle will warm slightly. After the initial charging period, you may begin using your system.


2. Install or Remove Attachment.

Select a cleaning attachment. Push straight onto the end of the power handle until you hear a "click". You may face the attachment position anyway that feels comfortable. DO NOT try to pull the attachment straight off. To remove, push or "bend" where the attachment meets the handle while pulling off at the same time.

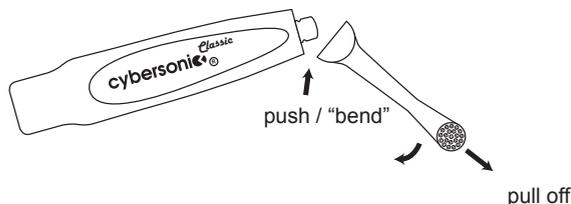

3. Power On.

Locate the push button on the side of the power handle. Push this button to power on and again to power off when you're finished cleaning. The power turns off automatically when the power handle is placed into the charging stand. The handle is designed so you can use it anyway that feels comfortable in your hand.

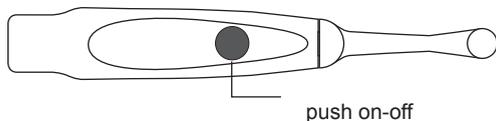

Charging the Batteries

Keep the charging stand plugged in continuously. For optimum performance you should place the power handle back into the charging stand after every use. You cannot overcharge the batteries. To extend battery life, unplug the charging stand every 2 months and plug the charger stand back in after 4 days or when the handle runs out of power. Batteries should last three to five years depending on how frequently the product is used.

Battery Disposal

We use nickel metal hydride rechargeable batteries. Please do not dispose of this appliance in the household wastebasket. Instead, dispose at a recycling center or according to local regulations.

Cleaning the System

After every use, rinse and clean any toothpaste film or residue from the attachments and the power handle. This is especially important if you use baking soda-peroxide toothpastes. Remove the cleaning attachments and rinse separately, then wipe them dry.




Replacing the Brush Attachments

Dentists recommend replacing your brush every 3 months. Replacement brushes may be purchased online at www.sonictoothbrush.com or by calling toll-free at 800-770-7525. All Cybersonic® customers are eligible for our FREE Brush Replacement Program. Just pay a small S&H fee, and we will ship you FREE brushes every 3 months for as long as you own your Cybersonic® Classic. Any active replacement program also includes a FREE Lifetime Warranty that covers the Charging Base and Power Handle. Now, Cybersonic® is truly the last toothbrush you'll ever have to buy!

Using the System

The System is an advanced sonic oral-care system that disintegrates plaque on contact and features a full range of attachments for complete oral health. Follow the simple instructions below to begin using the attachments properly.

Sonic Brushing Attachments

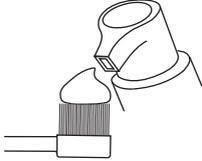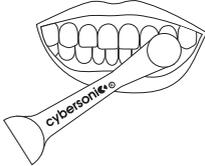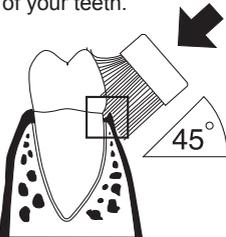- Place the Sonic Brush attachment onto the power handle. Remove the brush head cap and apply toothpaste or gel if desired. 
- Using moderately firm pressure, place the bristles on your teeth. Turn on the power handle. Move the brush head evenly over all surfaces of your teeth. No extra scrubbing or brushing motion is required in brushing your teeth and gums. Just press firmly and the sonic technology does the rest. 
- When brushing the sides of your teeth, place the brush head at a 45° angle. Apply firm pressure while moving along the gum line. The "Modified-Bass" technique directs sonic cleaning below the gum line. Spend at least two minute brushing all surfaces of your teeth. 


Use Moderate-Firm Pressure



Sonic Cleaning Energy Below The Gum Line

Tongue Cleaning Attachments

- Place the Sonic Tongue Cleaner attachment onto the power handle. You can use it with or without toothpaste or gel.

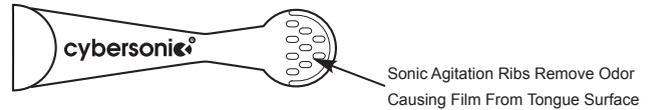

- Turn on power and place tongue cleaner on top of your tongue. Try to clean as far back towards the throat as possible. Pull or scrape the residue out and expel.

Sonic Flosser Attachment and Pre-Strung Replacement Floss Cartridges

- Place the flosser attachment onto the power handle. Then snap the pre-strung floss cartridge into the flosser attachment. 
- Turn on the power handle. Insert floss between teeth. For tight teeth, you may need to move the flosser from side to side.
- To clean between teeth, move flosser up and down both sides of your teeth. Also, floss under the gum line.
- To eject the floss cartridge, place thumbnail under edge then snap upward and away. 
- Dispose floss cartridge after each use. To order replacements go to www.sonictoothbrush.com or you can call customer service toll free at 1-800-770-7525.

Sonic White® One-Step Bleaching Gels (optional)

- Apply gel onto the Sonic Bleaching Brush head (optional). Push the plunger in until it reaches the one dosage marker and replace end cap on syringe. Each mark equals one day's treatment  6 Day Treatment (skip day 7, no treatment)
- Spread an even film of gel across FRONT teeth. Turn on power handle and make small circular motions. Continue for one to two minutes. Rinse. Look for results within one to three weeks.

Note: One syringe contains a one-week treatment; six days on, one day off. A one-week treatment is enough for normal yellowing or normal stains. Continue for additional 1-2 weeks for more stubborn stains. Sonic White is safe for daily use. For best results, do not bleach teeth right after using a fluoride toothpaste or rinse. During normal use, minor gum or tooth sensitivity may occur or white film or spots may appear, but are harmless and temporary. If you experience prolonged sensitivity, pain or bleeding gums, stop use and consult your dentist.

Commonly Asked Questions

- Q. Will the system loosen crowns, fillings, caps, bondings, veneers or other permanent dental work?
- A. No, Cybersonic® is designed to be used with all dental work.
- Q. How often should I replace my brush heads?
- A. Dentists recommend changing every three months. Your brush may not look worn but old brush bristles do not clean as well.
- Q. Is it safe for a person with a pacemaker?
- A. It is always advisable that a person with a special medical condition, including pacemakers, consult with a physician or pacemaker manufacturer.
- Q. Can I travel with my power toothbrush?
- A. Yes. A fully charged brush will hold a charge for at least one week of normal use, and the brush head cap keeps the bristles clean.
- Q. Is it suitable for kids?
- A. Yes! The Original round brush is actually a perfect size for children. You may start your kids using the system at the same age you would a manual toothbrush. Parental supervision is always advised with a new oral cleaning product.
- Q. Why do I sometimes need to push the ON/OFF button more than once?
- A. The electronic switch takes 1 second to re-set before it will respond to another push. Simply wait 1 second before pushing the button again.
- Q. How do I purchase replacement attachments?
- A. Replacement attachments may be purchased online at www.sonictoothbrush.com or by calling customer service toll-free at 800-770-7525. You may also sign up for our FREE Brush Replacement Program via the same toll free number or website. Pay a small S&H fee and we'll send you FREE replacement brushes directly to your home every 3 months.

Product contains nickel-metal hydride rechargeable batteries. Must be recycled or disposed of properly.



Important Safeguard Precautions!

Please Read Before Use.

Our toothbrushes have been carefully engineered for effectiveness and safety. However, basic safety precautions should always be followed when using any electrical product.

To reduce the risk of electrical shock, fire or injury:

- Avoid placing the charging stand where it might fall or be pulled into a tub, shower, sink, or toilet.

- Do not use the system while bathing or showering.
- Do not reach for charger if submerged in water. Unplug the stand immediately and drain the water before retrieving it by the cord.
- Keep the cord away from heated surfaces.
- Do not use the system if the charging stand, cord, or plug are damaged or have been immersed in water. Call Customer Service at 1-800-770-7525 for further instructions.
- Supervise children under 12 and disabled people when they are using this electrical device.

Warnings

- Avoid using toothpaste that contains peroxide and baking soda which can react with the plastic and cause cracks or breakage. If you do use peroxide baking soda toothpaste(s), you must rinse and clean more thoroughly after each usage.
- Consult your dentist before using product if you've had periodontal disease, dental abscesses, oral infections, oral pain, oral discomfort, bleeding gums, loose teeth or toothaches.
- Consult your physician before use if you are a cardiac patient with a pacemaker and have concerns. However, the system produces only a small magnetic field in the power handle, about the same as the magnetic field emitted by an electric shaver.
- If you experience persistent pain, bleeding or discomfort during or after use, consult your dentist.
- Not intended to replace regular dental check-ups or cleanings.